

## UNITED STATES AIR FORCE RESEARCH LABORATORY

---

### TESTING AND EVALUATION OF THE PROTOCOL SYSTEMS, INC., PROPAQ 206 EL ENCORE VITAL SIGNS PATIENT MONITOR

Edward W. Hade

KRUG Life Sciences, Inc.  
2504 Gillingham Drive, Suite 25  
Brooks AFB, Texas 78235-5104

James Sylvester, Major, USAF, NC

HUMAN EFFECTIVENESS DIRECTORATE  
FLIGHT STRESS PROTECTION DIVISION  
SYSTEMS RESEARCH BRANCH  
2504 Gillingham Drive, Suite 25  
Brooks AFB, Texas 78235-5104

July 1998

*Approved for public release; distribution is unlimited.*

19981218 045

## NOTICES

This final technical report was submitted by personnel of the Systems Research Branch, Crew Technology Division, Air Force Research Laboratory, AFMC, Brooks Air Force Base, Texas, under job order 7184-56-01.

This report was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor any agency thereof, nor any of their employees, nor any of their contractors, subcontractors, or their employees, makes any warranty, expressed or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency, contractor, or subcontractor thereof. The views and opinions of the authors expressed herein do not necessarily state or reflect those of the United States Government or any agency, contractor, or subcontractor thereof.


When Government drawings, specifications, or other data are used for any purpose other than in connection with a definitely Government-related procurement, the United States Government incurs no responsibility or any obligation whatsoever. The fact that the Government may have formulated or in any way supplied the said drawings, specifications, or other data, is not to be regarded by implication, or otherwise in any manner construed, as licensing the holder or any other person or corporation; or as conveying any rights or permission to manufacture, use or sell any patented invention that may in any way be related thereto.

The Office of Public Affairs has reviewed this report, and it is releasable to the National Technical Information Service, where it will be available to the general public, including foreign nationals.


This report has been reviewed and is approved for publication.

Government agencies and their contractors registered with Defense Technical Information Center (DTIC) should direct requests for copies to: Defense Technical Information Center, 8725 John J. Kingman Rd., STE 0944, Ft. Belvoir, VA 22060-6218.

Non-Government agencies may purchase copies of this report from: National Technical Information Services (NTIS), 5285 Port Royal Road, Springfield, VA 22161-2103.

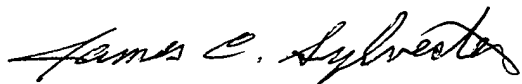

JAMES C. SYLVESTER, Major, USAF, NC  
Chief, A. F. Medical Equip & Dev. Laboratory



ROGER L. STORK, Colonel, USAF, BSC  
Chief, Flight Stress Protection Division

| REPORT DOCUMENTATION PAGE  |   |  | Form Approved<br>OMB No. 0704-0188   |  |
|--|---|--|--|--|
| Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503. |   |  |  |  |
| 1. AGENCY USE ONLY (Leave blank)   |   | 2. REPORT DATE<br>July 1998                                    |  | 3. REPORT TYPE AND DATES COVERED<br>Final, July 1998 |
| 4. TITLE AND SUBTITLE<br>Testing and Evaluation of the Protocol Systems, Inc. PROPAQ 206 EL Enclore Vital Signs Patient Monitor  |   |  | 5. FUNDING NUMBERS<br>PE: 62202F<br>PR: 7184<br>TA: 56<br>WU: 01           |  |
| 6. AUTHOR(S)<br>Edward W. Hade<br>James Sylvester, Major, USAF   |   |  |  |  |
| 7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES)<br>Air Force Research Laboratory (AFMC)<br>Human Effectiveness Directorate<br>Flight Stress Protection Division<br>2504 Gillingham Drive, Suite 25<br>Brooks AFB, TX 78235-5104   |   |  | 8. PERFORMING ORGANIZATION<br>REPORT NUMBER<br><br>AFRL-HE-BR-TR-1998-0053 |  |
| 9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)  |   |  | 10. SPONSORING/MONITORING<br>AGENCY REPORT NUMBER                          |  |
| 11. SUPPLEMENTARY NOTES  |   |  |  |  |
| 12a. DISTRIBUTION AVAILABILITY STATEMENT<br><br>Approved for public release; distribution is unlimited.  |   |  | 12b. DISTRIBUTION CODE   |  |
| 13. ABSTRACT (Maximum 200 words)<br><br>The Protocol Systems Inc., 206 EL Monitor is a portable, self-contained, general purpose monitor capable of continuous reception and display of multiple patient physiological parameters. The unit operates on 115 VAC/60-400 Hz, external 12 VDC, and an internal battery pack. The unit weighs approximately 12.68 lbs and its dimensions with the expansion module are 9.65 in. H X 8.25 in. W X 7.56 in. D. It has a separate UPA/Style B 503-0054-00 power adapter as an available attachment.   |   |  |  |  |
| 14. SUBJECT TERMS<br>Aircraft      Aeromedical      Protocol<br>206 EL      Medical equipment<br>Airworthy      physiological  |   |  | 15. NUMBER OF PAGES<br>26  |  |
|  |   |  | 16. PRICE CODE   |  |
| 17. SECURITY CLASSIFICATION<br>OF REPORT<br><br>Unclassified   | 18. SECURITY CLASSIFICATION<br>OF THIS PAGE<br><br>Unclassified | 19. SECURITY CLASSIFICATION<br>OF ABSTRACT<br><br>Unclassified | 20. LIMITATION OF<br>ABSTRACT<br><br>UL                                    |  |

**DTIC QUALITY INSPECTED 3**

## TABLE OF CONTENTS

|  |    |
|--|----|
| BACKGROUND .....                               | 1  |
| DESCRIPTION.....                               | 1  |
| PROCEDURES .....                               | 2  |
| INITIAL INSPECTION AND TEST PREPARATION.....   | 3  |
| TEST SETUP .....                               | 3  |
| PERFORMANCE CHECK .....                        | 5  |
| VIBRATION.....                                 | 5  |
| ELECTROMAGNETIC COMPATIBILITY .....            | 7  |
| THERMAL/HUMIDITY ENVIRONMENTAL CONDITIONS..... | 8  |
| HYPOBARIC CONDITIONS .....                     | 9  |
| AIRBORNE PERFORMANCE .....                     | 10 |
| EVALUATION RESULTS .....                       | 10 |
| INITIAL INSPECTION.....                        | 10 |
| VIBRATION.....                                 | 10 |
| ELECTROMAGNETIC COMPATIBILITY .....            | 10 |
| THERMAL/HUMIDITY ENVIRONMENTAL CONDITIONS..... | 11 |
| HYPOBARIC CONDITIONS .....                     | 11 |
| AIRBORNE PERFORMANCE .....                     | 11 |
| SUMMARY.....                                   | 11 |
| REFERENCES .....                               | 13 |
| APPENDIX.....                                  | 14 |

## LIST OF FIGURES

|  |   |
|--|---|
| Figure 1. The Protocol Systems, Inc., Propaq 206 EL Encore .....         | 2 |
| Figure 2. Test Setup .....   | 4 |
| Figure 3. Vibration Table Mounting.....                                  | 5 |
| Figure 4. MIL-STD-810E, Category 10, figures 514.4-16 and 514.4-17 ..... | 6 |

## **ACKNOWLEDGMENTS**

We would like to thank those who helped and provided advice during the evaluation of the Propaq 206 EL Encore. We would particularly like to thank:

MSgt Butch Blake  
TSgt Allen Jones  
Mr. Douglas Townsend

**TESTING AND EVALUATION OF THE  
PROTOCOL SYSTEMS, INC.  
PROPAQ 206 EL ENCORE  
VITAL SIGNS PATIENT MONITOR**

**BACKGROUND**

Representatives of Protocol Systems, Inc. requested an Aeromedical Research evaluation of their Propaq 206 EL Encore vital signs monitor for use on board USAF aeromedical evacuation aircraft.

**DESCRIPTION**

The unit tested was the Protocol Propaq Encore model 206 EL, SN: DA006428 with expansion module SN: DC003733. The Propaq Encore expansion module attaches to the monitor and houses additional capabilities. This expansion module was fitted with Printer, Pulse Oximetry (SpO<sub>2</sub>), and Capnography (CO<sub>2</sub>) options. The unit tested will hereby be referred to as the 206 EL, or Encore (figure 1). This unit is a light-weight portable patient monitor capable of monitoring: ECG (1 channel: 3-lead\*); NIBP, noninvasive blood pressure, (1 channel: cuff); IBP, invasive blood pressure, (2 channels); temperature (2 channels: YSI-400 and 700 series-compatible connectors); pulse oximetry (1 channel: SpO<sub>2</sub>); CO<sub>2</sub> (1 channel); and respiratory rate. This unit has a printer and Hewlett Packard Connector-Compatible Side Panel. The display in the 206 EL is Electroluminescent (EL). With printer/SpO<sub>2</sub>/CO<sub>2</sub>, the dimensions of the unit are as follows: height, 9.65 in. (24.5 cm); width, 8.25 in (20.9 cm); depth, 7.56 in (19.2 cm); weight, 12.68 lb (5.8 Kg). The 206 EL has an internal, 8 V/6 amp-hr, sealed gel-type lead-acid battery. Battery life is rated at 3.5-4.5 hours depending on product configuration with a recharge time of 8-12 hours with the instrument on, or 6 to 8 hours with the instrument off. The unit has an adapter, which converts 100-120 VAC/60Hz to 16-24 VDC/25 VA: part number 503-0054-00. The unit can also be powered by an external 12-28 VDC source.

\*Note: This device also offers the option of a 5 lead ECG monitoring system. This 5 lead option was not evaluated as part of airworthiness testing.

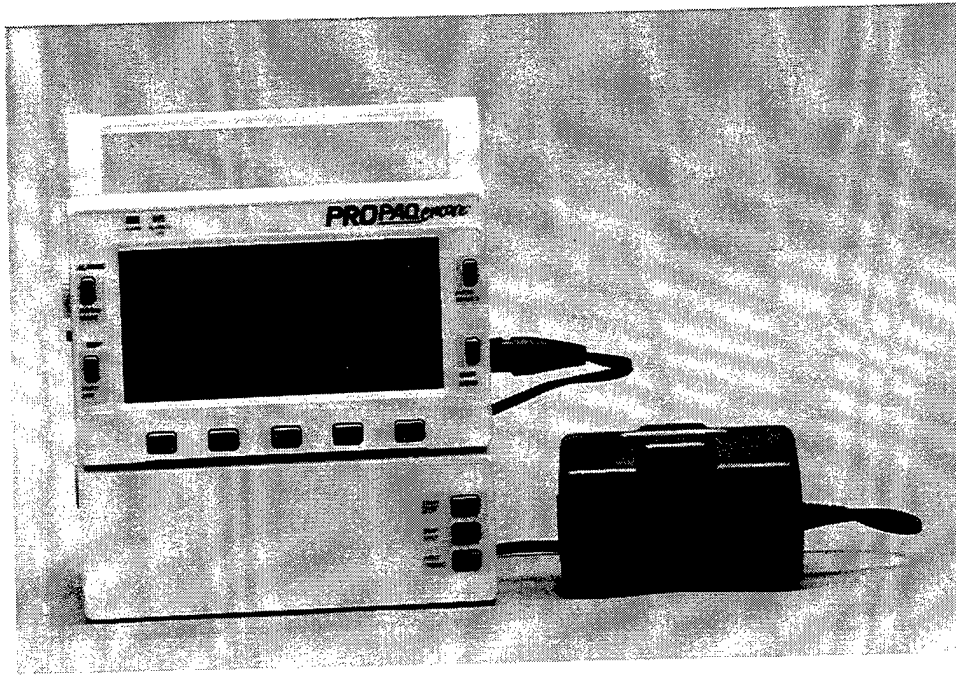


Figure 1. The Protocol Systems, Inc., Propaq 206 EL Encore.

### PROCEDURES

Test methods and performance criteria were derived from various military standards (Reference List 1-4), nationally recognized performance guidelines (5), and manufacturer's literature (6). The Aeromedical Research Procedures Guide describes additional safety and human interface issues to be considered during equipment testing (7). A test setup and performance check were developed specific to this product to verify proper functioning of the equipment when subjected to various tests representing the airborne environment and stresses of flight.

The device was subjected to the following laboratory and inflight tests to observe and evaluate its performance under anticipated operational conditions.

1. Initial Inspection
2. Vibration
3. Electromagnetic Compatibility (EMC)

4. Thermal/ Humidity Environmental Conditions, encompassing:

- a. Hot Operation
- b. Cold Operation
- c. Humidity
- d. Hot Temperature Storage
- e. Cold Temperature Storage

5. Hypobaric Conditions

- a. Cabin Pressure/Altitude
- b. Rapid Decompression to Ambient pressure

6. Airborne Feasibility

**INITIAL INSPECTION AND TEST PREPARATION**

a. The Propaq 206 EL Encore was inspected for quality of workmanship, production techniques and possible damage that might have occurred during shipment.

b. The Propaq 206 EL Encore was checked to ensure it met safety requirements and operating characteristics established in National Fire Protection Agency (NFPA) 99, Standards for Health Care Facilities (7), Electrical Shock Hazards, AFI 41-203 (8), and Equipment Management in Hospitals, AFI 41-201 (9). Ground resistance and leakage current measurements were made at 115 VAC/60 Hz.

c. The Propaq 206 EL Encore was examined to verify it met basic requirements for acceptable human factors design as outlined in MIL-STD 1472 (3).

d. A test setup and performance check were developed to evaluate the Propaq 206 EL Encore's operation in accordance with manufacturer/customer specifications throughout the various testing conditions.

**TEST SETUP**

a. Connect the three ECG leads from the Encore to the corresponding (color coded) receptacles on the Lionheart.

- b. Plug the YSI temperature cable (1/4 inch phone jack side) into the Encore's "T1" port. Plug the opposite end of the cable (tri-axil side) into the Lionhearts' 700 series temperature output connector.



c. Configure the Lionheart with the following settings:

- Temperature: 30°C
- Lead Select: III
- ECG Amplitude: 1.0 mV

d. Secure the non-invasive tubing line to the NIBP port on the Encore. Using a T connector, attach the Cufflink inline between the BP cuff and the non-invasive tubing attached to the NIBP port. Wrap the BP cuff tightly around the appropriate adult cuff mandrel. For an adult cuff, use two end and two spacer blocks. After zeroing the transducer, the Cufflink is configured to: ADAMS Adult-120/80 (90).

e. Plug the SpO<sub>2</sub> cable into its corresponding port on the Encore (9 pin connector) and the other end into the Nellcor Pulse Oximeter Simulator. Set the Nellcor Pulse Oximeter Simulator to 98% SpO<sub>2</sub> and a pulse rate of 60 bpm.

f. Plug the CO<sub>2</sub> sensor cord into its corresponding port on the Encore. Attach the sensor to the CO<sub>2</sub> line. Secure it to one end of a section of corrugated ventilator tubing. At the other end, place a moisture trap filter.

g. To limit paper usage during evaluation, loop a four inch piece of printer paper through the printer and secure the ends with tape.

The Encore will continually monitor temperature, SpO<sub>2</sub>, and CO<sub>2</sub>. The NIBP operation can be initiated manually or programmed at set intervals.

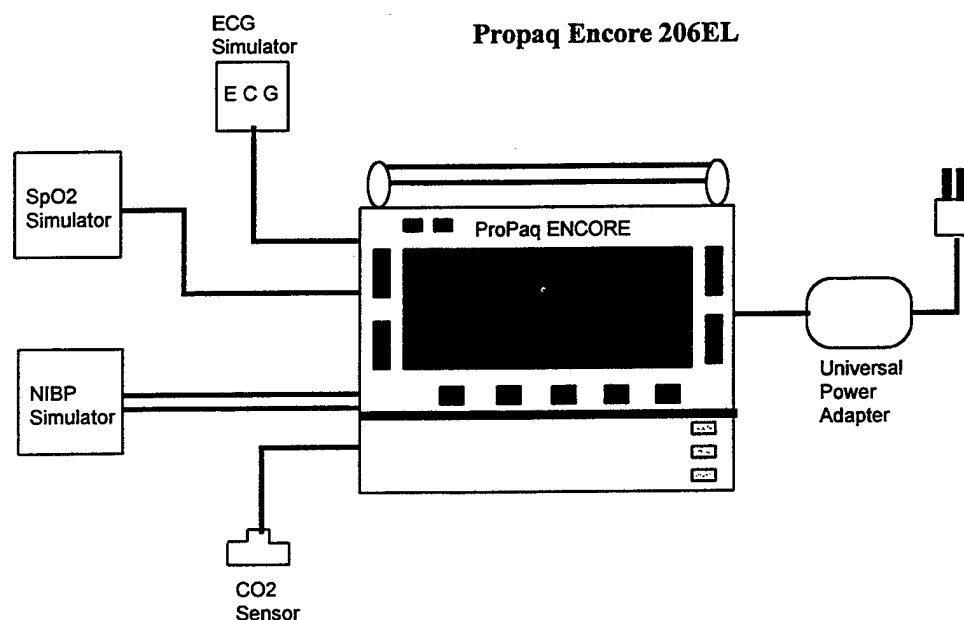


Figure 2. Test Setup

## **PERFORMANCE CHECK**

The Performance Check, as outlined in the approved test plan, was used to validate the function of the 206 EL in each of the test conditions. Measurements were taken during initial operation at standard ambient conditions and served as a baseline for later comparison. The performance check consisted of recording the values for each monitored physiologic parameter three times and activating the printer to ensure its function. In many cases, the 206 EL was continuously monitored through the duration of the test. Performance checks occurred at defined intervals throughout the test.

## **VIBRATION**

Vibration testing is critical to determine "the resistance of equipment to vibrational stresses expected in its shipment and application environments" (2). Vibration testing was conducted at Aeromedical Research's vibration facility. This testing involved a set of operational tests performed along each of the Encore's three axes - X, Y, and Z; the Encore's components were mounted on the NATO litter segment on the vibration table as they would be in the aircraft (Figure 3). They were subjected to vibration curves with levels and lengths derived from MIL-STD-810E, Category 10, Figures 514.4-16 and 514.4-17 (Figure 4).

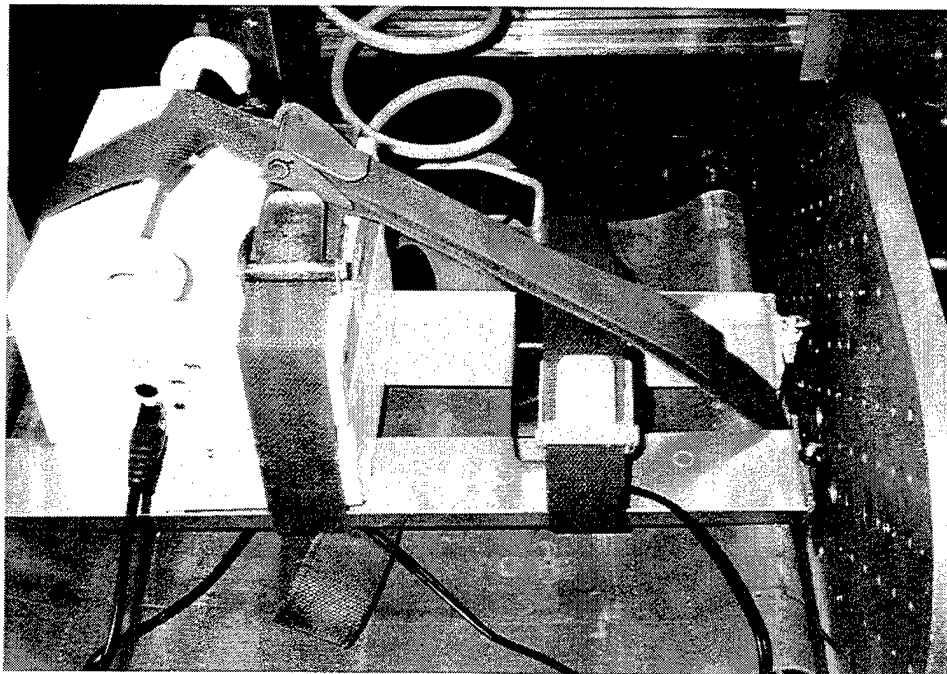


Figure 3. Vibration Table Mounting

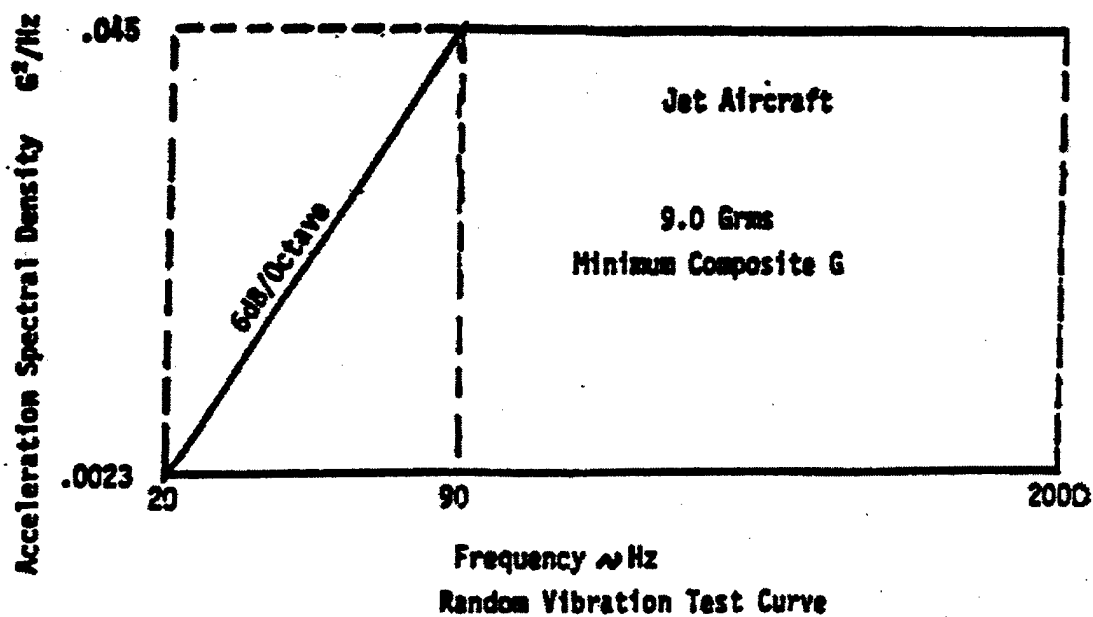
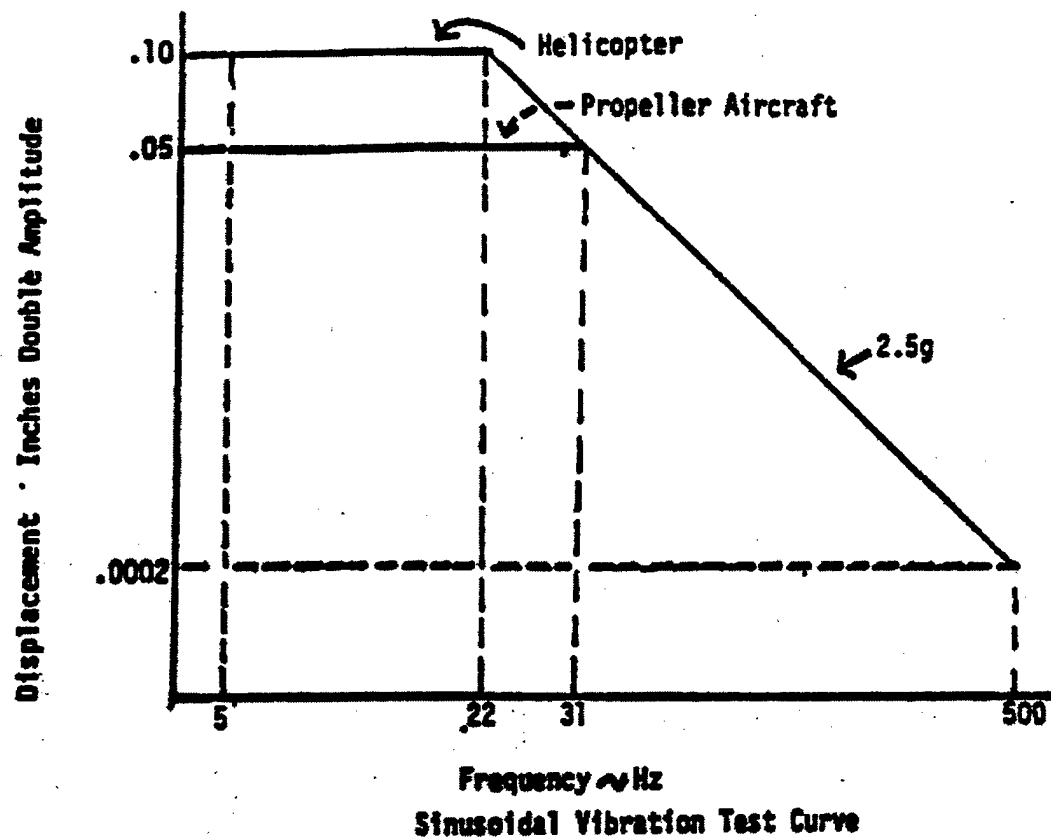


Figure 4. MIL-STD-810E, Category 10, figures 514.4-16 and 514.4-17

## **ELECTROMAGNETIC COMPATIBILITY**

Electromagnetic compatibility testing is a primary concern on USAF aeromedical evacuation aircraft. Ensuring the safety of everyone on board is the driving factor to assessing the effects of excessive electromagnetic emissions and their influence on aircraft navigation and communication equipment. Medical devices may also be susceptible to fields generated by the aircraft equipment or other medical devices and malfunction in their presence.

The Encore was evaluated for compliance with MIL-STD-461D (1) and 462D (2) . ASC/ENAI, Wright-Patterson AFB performed all of the EMI evaluation in their electromagnetic compatibility facility and determined the airworthiness of the medical device. Specific tests conducted were as follows:

a. Radiated Emissions (RE-102), "Radiated Emissions, Electric Field, 10 kHz to 18 GHz.": For Air Force aircraft applications, radiated emissions were tested in a narrower range of frequencies from 2 MHz - 1 GHz. This test measured the amount of EMI emitted by the equipment during its operation. It was performed to verify that the device does not affect other pieces of equipment that may be susceptible to electromagnetic emissions (i.e., aircraft navigation and communication equipment).

b. Conducted Emissions (CE-102), "Conducted Emissions, Power Leads, 10 kHz to 10 MHz.": For Air Force aircraft applications, conducted emissions were tested throughout the entire band of 10 kHz - 10 MHz. This test which measured emissions generated by the medical device along its power supply lines, was performed to verify that operating the device using line power does not affect other items connected to the same power source, particularly aircraft systems.

c. Radiated Susceptibility (RS-103), "Radiated Susceptibility, Electric Field, 10 kHz to 40 GHz.": For Air Force aircraft applications, radiated susceptibility was tested in a narrower frequency range from 30 MHz - 12.4 GHz at the following field strength levels: 20 V/M below 1 GHz and 60 V/M above 1 GHz (field strength values from Table IV, category Aircraft Internal, of 461D). This test determined whether or not the device would withstand pre-defined levels of EMI generated by antennas both internal and external to the aircraft.

d. Conducted Susceptibility (CS-101), "Conducted Susceptibility, Power Leads, 30 Hz to 50 kHz.": For Air Force aeromedical aircraft applications, conducted susceptibility was tested throughout the entire frequency band, from 30 Hz to 50 kHz. This test determined whether the components would "withstand ripple voltages associated with allowable distortion of power source voltage wave forms."

e. Conducted Susceptibility (CS-114), "Conducted Susceptibility, Bulk Cable Injection, 10 kHz to 400 MHz.": For Air Force aeromedical aircraft applications conducted susceptibility was tested throughout a narrower portion of the frequency band, from 10 kHz to 200 MHz. This test determined whether "simulated currents that will be developed on platform cabling from electromagnetic fields generated by antenna transmission would affect the equipment under test."

f. Conducted Susceptibility (CS-115), "Conducted Susceptibility, Bulk Cable Injection, Impulse Excitation": This test was performed to verify that the Encore could withstand the "fast rise and fall time that may be present due to platform switching operations and external transient environments such as lightning and electromagnetic pulse."

g. Conducted Susceptibility (CS-116), "Conducted Susceptibility, Damped Sinusoidal Transients, Cables and Power leads, 10 kHz to 100 MHz," respectively. The "basic concept of this test is to simulate electrical current and voltage waveforms occurring in platforms from excitation of natural resonances."

During emissions testing, all options were operating for the duration of the test to create the "worst case" emissions scenario. Throughout the testing, the recorder (printer) ran continuously, and the apnea alarm continuously sounded at maximum volume. The 206 EL was in turbo-cuff mode, such that the NIBP option was continuously activated. For susceptibility testing, the unit was operated as described earlier in the equipment set-up and performance check sections. For both emissions and susceptibility testing, the 206 EL was tested for operation on 115 VAC/60 Hz, and internal batteries.

### **THERMAL/HUMIDITY ENVIRONMENTAL CONDITIONS**

Extreme temperature and humidity testing is critical to determine if aeromedical equipment can be stored and operated under severe environmental conditions "without experiencing physical damage or deterioration in performance" (3). Extreme environmental conditions can have numerous detrimental effects on medical equipment including, but not limited to: changes in material characteristics and material dimensions, possible overheating, changes in lubricant viscosity, changes in electronic components, and electronic or mechanical failures due to rapid water or frost formation.

Testing was conducted in the Armstrong Laboratory's, Thermotron Industries, Model SM-32C environmental chamber operated and monitored by aeromedical research personnel assigned to the Systems Research Branch (CFTS) of the Crew Technology Division at Armstrong Laboratory, Brooks AFB, TX. The 206 EL was placed in the center of the environmental chamber. All input and output cables and wires were routed through a port in the chamber wall, which was subsequently sealed with a precut sponge plug. The other components of the test setup were outside the chamber. For operational tests, the 206 EL was monitored continuously and a performance check was conducted every fifteen minutes. For storage tests, the 206 EL was placed in the chamber and remained nonoperational throughout the storage portion of the test. Upon completion of this test the chamber and device was brought to standard ambient conditions for 30 minutes. Aeromedical Research personnel then conducted a performance test and monitored the unit for one hour to verify successful operation. The following describe the conditions of the environmental tests performed:

- a. Humidity:  $94 \pm 4\%$  RH,  $85^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$  ( $29.5^{\circ}\text{C} \pm 2^{\circ}\text{C}$ ) for 4 hr
- b. Hot Temp Operation:  $120^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$  ( $49^{\circ}\text{C} \pm 2^{\circ}\text{C}$ ) for 2 hr
- c. Cold Temp Operation:  $32^{\circ}\text{F} \pm 7.2^{\circ}\text{F}$  ( $0^{\circ}\text{C} \pm 4^{\circ}\text{C}$ ) for 2 hr
- d. Hot Temp Storage:  $140^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$  ( $60^{\circ}\text{C} \pm 2^{\circ}\text{C}$ ) for 6 hr
- e. Cold Temp Storage:  $-40^{\circ}\text{F} \pm 3.6^{\circ}\text{F}$  ( $-40^{\circ}\text{C} \pm 2^{\circ}\text{C}$ ) for 6 hr

### **HYPOBARIC CONDITIONS**

Testing was conducted in the Armstrong Laboratory research chambers which were operated and monitored by chamber operations personnel assigned to the Systems Research Branch (CFTS) of the Crew Technology Division at Armstrong Laboratory, Brooks AFB, TX.

a. Hypobaric Chamber Testing: Altitude testing is critical for aeromedical evacuation equipment due to the potential effects of barometric pressure changes on the equipment. A majority of the aircraft that are characterized as opportune aircraft available for use in aeromedical evacuation pressurize their cabin to barometric pressures equivalent to 8,000-10,000 feet above sea level. However, the differences in pressures can be critical to the effective operation of some medical equipment. Altitude testing consisted of operating the 206 EL while ascending from ground level to 10,000 feet (maintaining altitude for one hour) and then descending back to ground, at rates of 5000 ft/min, while stopping at 2000 ft increments to allow for performance checks.

b. Rapid Decompression Testing: Rapid decompressions are caused by the loss of aircraft cabin pressurization and subsequent pressure equalization with ambient atmospheric pressures. It is important to determine how medical equipment will function during and after such a decompression and verify that it will not endanger a patient, the aircraft personnel, or the aircraft itself. The 206 EL operated inside the rapid decompression test chamber as the chamber was depressurized to an equivalent of 8,000 ft altitude. Then, the chamber altitude was brought to 45,000 ft over a period of 60 seconds, held at 45,000 ft for a few minutes, and then brought back down to ground at a rate of 10,000-12,000 ft/min. The test was repeated twice with the decompressions occurring over seven and one seconds, respectively. The 206 EL was monitored throughout the series of decompressions, including performance checks each time the unit returned to ground level. The simulator equipment remained outside the chamber. Cables joining the Lionheart, Cufflink, and Nellcor to the 206 EL were run through putty-sealed access ports in the chamber walls.

## **AIRBORNE PERFORMANCE**

Airborne performance evaluations are a cost-effective and invaluable means of validating equipment clinical and operational performance during actual operating conditions. By carefully evaluating medical equipment items in their actual environment, Aeromedical Research ensures that all pertinent patient care issues are adequately addressed by the test protocols. Safe and reliable operation of this medical device is the primary goal of the inflight evaluation and forms the basis for subsequent recommendations to the users.

This phase of testing was conducted by an aircraft-qualified aeromedical flight nurse and aeromedical research technicians on board both a C-9 and C-130 aeromedical evacuation mission. The 206 EL was secured to the litter and evaluated throughout the flights by Aeromedical Research technicians as well as by the other members of the aeromedical evacuation crew. Human factors characteristics, securing methods, and equipment setup times and locations were also evaluated.

## **EVALUATION RESULTS**

### **INITIAL INSPECTION**

Initial inspection revealed no manufacturing defects. The unit performed to the manufacturer's specification. Electrical safety test results showed that all parameters were within referenced guideline limits.

### **VIBRATION**

The Propaq Encore operated within manufacturer's specifications throughout vibration testing.

### **ELECTROMAGNETIC COMPATIBILITY**

While at the Wright-Patterson AFB EMI testing facilities, the Propaq Encore experienced EMI failures. After initial electromagnetic interference evaluations, ASC/ENAE, Wright Patterson AFB, approved the Propaq Encore model 206 EL for use on large-bodied USAF aircraft only. The Wright-Patterson EMI certification letter stated that the 206 EL would not interfere with communication or navigation systems on large bodied aircraft, and it would be unlikely that the normal aircraft RF transmitting systems would cause interference to the 206 EL with the following exceptions: "heart rate" and "ECG" may be affected and should not be solely relied upon for clinical judgements in critical situations. Protocol Systems, working with Wright-Patterson AFB, WL/AASW, and Aeromedical Research, Brooks AFB developed new circuitry within the 206 EL (Printer and CO2) in order to overcome these EMI difficulties.

WL/AASW tested and ASC/ENAE certified the newly modified 206 EL in June 1997 for operation during all phases of flight on all USAF aircraft.

### **THERMAL/HUMIDITY ENVIRONMENTAL CONDITIONS**

The Propaq Encore model 206 EL operated satisfactorily during all five phases of testing. Testing was conducted in the Armstrong Laboratory's Thermotron Environmental Chamber operated and monitored by Aeromedical Research personnel assigned to the Systems Research Branch (CFTS) of the Crew Technology Division at Armstrong Laboratory, Brooks AFB, TX.

### **HYPOBARIC CONDITIONS**

1. Cabin Pressure/Altitude: The Propaq Encore model 206 EL performed in accordance with manufacturer's specifications throughout testing.
2. Rapid Decompression: The 206 EL operated within manufacturer's specifications following each rapid decompression and did not present a safety hazard throughout the decompression.

### **AIRBORNE PERFORMANCE**

The inflight evaluation of the Propaq Encore model 206 EL was performed on a C-9 aeromedical evacuation mission and C-130 aeromedical readiness mission. The inflight evaluations of the 206 EL were successfully completed with the following comments: (1) the audible alarms are difficult to hear in the noisy aircraft environment, and (2), the alarm indications are difficult to view from the side of the unit. For these reasons, Aeromedical Research recommends that the 206 EL be mounted such that a crew member is monitoring the display from a front view. The securing capabilities with the 206 EL are adequate, utilizing litter equipment brackets with litter straps or the Waters bracket. However, patient connectors and the display limit positioning of the cargo straps. Analysis of flight data indicated that this unit was easy to enplane and deplane and was compatible with aircraft electrical systems.

### **SUMMARY**

The test and evaluation of Protocol System's Propaq 206 EL, SN: DA006428, and expansion module, SN: DC003733 is complete. Aeromedical Research found this unit, and all 206 ELs with serial numbers higher than EA000225\*, acceptable for use during all phases of flight on all USAF aircraft (including small and large body, fixed and rotary wing) while operating on battery, 115 VAC/60 Hz, and 28 VDC in the aeromedical evacuation environment with the following comments and recommendations:



a. Because the carbon dioxide and breath rate sensor ceased operation during the laboratory's hot operation test (120°F), Aeromedical Research recommends restricting operational use in extreme hot environments if the carbon dioxide and breath rate is a critical portion of patient monitoring.

b. Aeromedical Research recommends that the 206 EL be mounted such that a crew member is consistently able to monitor the display from a front view because the audible alarms are difficult to hear in the noisy aircraft environment and the visual alarm indicators are difficult to view from the side of the unit. The securing capabilities with the 206 EL proved adequate, utilizing litter equipment brackets with litter straps or the Waters bracket. However, patient connectors and the display limit positioning of the cargo straps may indicate the need for a more versatile mounting system for the 206 EL.

c. The 206 EL has many additional features/options to include: HP connectors, multiple power adapters, and defibrillator synchronization. The HP connector-compatible option makes the Propaq Encore compatible with many Hewlett-Packard sensors and accessories used with the Hewlett-Packard Component Monitoring System. This option was not tested and therefore not approved for use. As the Defibrillator Synchronization feature is designed to operate only with the Physio-Control LifePak 5 and LifePak 6 defibrillators, it was not tested and likewise not aeromedical airworthiness certified. The Propaq Encore is approved for use only with the UPA/Style B 503-0054-00 Power Adapter.

d. Based on the prior analysis of Protocol Systems comparison of the Propaq 206 EL with the 202 EL and 204 EL, Aeromedical Research has concluded that these physiologic patient monitors (with serial numbers higher than EA000225\*) will not require additional testing and can also be considered approved for use on all USAF aircraft.

e. Protocol Systems agreed to place a label on each Propaq Encore stating "Serial No. EA000225 and greater are approved for use during all phases of flight aboard U.S. Air Force aircraft." Label will be antique gold with black lettering. The label will be 3.69" W X 1.25" H X 0.25" Radius.

\* NOTE: Units with serial numbers lower than EA000225 were certified for operation during all phases of flight only on cargo (large body) USAF aircraft.

## REFERENCES

1. Department of Defense: *Requirements for the Control of Electromagnetic Interference Emissions and Susceptibility*. MIL-STD 461D, Washington DC: January 1993.
2. Department of Defense: *Measurement of EMI Characteristics*. MIL-STD 462D, Washington DC: February 1996
3. Department of Defense: *Environmental Test Methods and Engineering Guidelines*. MIL-STD 810E, Washington DC: 1989
4. Department of Defense: *Human Engineering Design Criteria for Military Systems, Equipment, and Facilities*. MIL-STD-1472D, Washington DC: March 1989
5. *Emergency Care Research Institute* (ECRI)
6. Protocol Systems, Inc., Propaq Encore 206 EL, Operator's Manuals.
7. *Aeromedical Research Procedures Guide*, Internal Operating Instruction, Systems Research Branch, Armstrong Laboratory.
8. National Fire Protection Agency (NFPA) 99, *Health Care Facilities Code*
9. AFI 41-203, *Electrical Shock Hazards*
10. AFI 41-201, *Equipment Management in Hospitals*

## APPENDIX

### MANUFACTURERS SPECIFICATIONS OF THE PROTOCOL SYSTEMS, INC. PROPAQ ENCORE 206 EL PHYSIOLOGIC PATIENT MONITOR

#### ECG SPECIFICATIONS

|                                |  |
|--------------------------------|--|
| CONNECTOR                      | AAMI 6 pin or Hewlett-Packard compatible 12-pin style connector (optional).  |
| SELECTABLE LEADS               | I, II, III, aVR, aVL, aVF, V   |
| LEAD FAULT INDICATOR           | LA, LL, RA, MULTIPLE   |
| ECG SIZE IN mV/cm              | 4, 2, 1, .5, .2  |
| DISPLAY SWEEP SPEEDS           | 12.5, 25, and 50 mm/sec  |
| QRS TONE VOLUME                | High, Low, Medium, Off   |
| QRS TONE FREQUENCY             | 900 Hz; variable pitch with SpO <sub>2</sub> option and SpO <sub>2</sub> being monitored   |
| FREEZE BUFFER                  | 3.9 sec at 25 mm/sec   |
| BANDWIDTH                      | 0.5 to 40 Hz   |
| INPUT PROTECTION               | Electrosurgery and defibrillator protected. All models also include electrosurgery interference suppression  |
| LEAD FAIL SENSE CURRENT        | 50 nA dc for active leads<br>100-200 nA dc for driven leads  |
| TALL T-WAVE REJECTION          | Meets and exceeds AAMI (USA) EC-1983, section 3.1.2.1, part 3, for 1.2 mV T-wave and 1 mV QRS using AAMI test waveform   |
| COMMON MODE REJECTION          | <1 mV p-p RTI for 10Vrms, 50/60 Hz input, input unbalanced, FILTER function OFF<br><.1 mV p-p RTI for 10 Vrms, 50/60 Hz input, input unbalanced, FILTER function ON                |
| INPUT IMPEDANCE                | >2.5 M differential at 60Hz  |
| INPUT RANGE (ac)               | +/- 5 mV   |
| INPUT RANGE (dc)               | up to +/-300mV   |
| QRS DETECTOR                   | Width range: 25-120 ms amplitude<br>Range: .3 to 5mV (RTI)   |
| HEART RATE COUNTER RANGE       | 25-250 bpm   |
| HEART RATE METER RESPONSE TIME | Responds to change in heart rate within 5-9 seconds depending on physiological waveform. (Including AAMI 3.1.2.1 parts 6 and 7 waveforms.) Includes 1-sec readout update interval. |
| HEART RATE ACCURACY            | +/- 3 bpm or 3%, whichever is greater  |
| HEART RATE AVERAGING METHOD    | see User's Guide   |
| DRIFT TOLERANCE                | 80 bpm indicated for 80 bpm ECG plus drift waveform  |
| PACER DISPLAY                  | Pacer indicator shown on screen if PACER function turned on; pacer spike always shown if of sufficient amplitude.  |

PACER PULSE REJECTION  
RESPONSES TO IRREGULAR RHYTHM

Ventricular Bigeminy (VB)  
Slowing Alternating VB  
Rapid Alternating VB  
Bidirectional Systole

see User's Guide

77-82 bpm  
63-81 bpm  
115-123 bpm  
87-93 bpm

## INVASIVE PRESSURE SPECIFICATIONS

TRANSDUCER TYPE  
TRANSDUCER EXCITATION IMPEDANCE  
RANGE

TRANSDUCER SENSITIVITY  
EXCITATION VOLTAGE  
CONNECTOR

BANDWIDTH  
ZERO DRIFT  
ZERO ADJUSTMENT  
NUMERIC ACCURACY

PRESSURE RANGE  
PULSE RANGE  
LEAKAGE CURRENT  
ELECTROSURGERY SUPPRESSION

Strain-gauge resistive bridge

200 - 2000  
5 micro V/V/mmHg  
4.85-V pulsed dc @ 181 Hz  
ITT-Cannon plug MS3106F-14S-6P Std.  
Hewlett-Packard compatible 12 in. connector  
digital filtered, dc to 25 Hz  
+/- 1 mmHg without transducer drift  
+/- 200 mmHg including transducer offset  
+/- 2mmHg or 2% of reading, whichever is  
greater, plus the transducer error  
-30 to 300 mmHg  
25-250 bpm  
Meets ANSI/AAMI risk requirements  
Included in all EL display monitors

## NIBP SPECIFICATIONS

METHOD  
CONTROL  
AUTO INTERVALS  
TURBOCUFF

DISPLAYED PRESSURES

SYSTOLIC RANGE  
DIASTOLIC RANGE  
MEAN RANGE  
NUMERIC ACCURACY  
MINIMUM INFLATION PRESSURE  
DEFAULT INFLATION PRESSURE  
CUFF OVERPRESSURE  
PULSE RATE RANGE

MAXIMUM DETERMINATION TIME  
TYPICAL DETERMINATION TIME  
TYPICAL DETERMINATION TIME WITH  
ARTIFACT  
MINIMUM TIME BETWEEN MEASUREMENTS  
ELECTROSURGERY SUPPRESSION

Oscillometric  
Automatic and manual measurement control  
1, 2, 3, 5, 10, 15, 30, and 60 min  
Maximum measurements allowable in a 5 min  
period  
Systolic, Diastolic, and Mean plus on-screen  
monitor  
30-260 mmHg  
20-235 mmHg  
25-255 mmHg  
+/- 3 mmHg  
100 mmHg  
Adult - 160 mmHg, Child - 120 mmHg  
280 mmHg  
30-220 bpm (without ECG)  
25-200 bpm (with ECG)  
3 min  
30-45 sec  
up to 70 sec  
30 sec  
Included in all EL display monitors

## PULSE OXIMETRY

|  |  |
|--|--|
| RANGE                                    | 0-100%   |
| PROBE ACCURACY<br>(specified at 28-42°C) | 70-100% +/- 2 digits, 0-70% unspecified                                      |
| PULSE RATE RANGE                         | 20-250 bpm   |
| PULSE RATE ACCURACY                      | +/- 3 bpm  |
| SENSOR COMPATIBILITY                     | Compatible only with NELLCOR sensors listed in Chapter 2 of the User's Guide |
| ELECTROSURGERY SUPPRESSION               | Included in all models 202 EL, 204 EL, 206 EL                                |

## CO<sub>2</sub> OPTION

|                        |   |
|------------------------|---|
| CO <sub>2</sub> SENSOR |   |
| Sensor type            | Mainstream  |
| Principle of operation | NDIR single-beam, single path/wavelength, ratiometric |
| Warm-up time           | 45 sec typical, 3 min maximum                         |
| Response time          | 30 ms typical, 60 ms maximum                          |
| Calibration            | Verify semi-annually, calibrate only as required      |

### CO<sub>2</sub> SENSOR AND CABLE DIMENSIONS AND WEIGHT

|               |                      |
|---------------|----------------------|
| Sensor Height | 1.003 in             |
| Sensor Width  | 1.036 in             |
| Sensor Depth  | .78 in               |
| Sensor Weight | < .39 oz             |
| Sensor Volume | 0.81 in <sup>3</sup> |
| Cable Length  | 10 ft nominal        |

### CO<sub>2</sub> AIRWAY ADAPTER

|           |  |
|-----------|--|
| Type      | Per ISO 3040, single-use                   |
| Size      | 15 mm ID (meets ISO specifications)        |
| Material  | clear polycarbonate, with sapphire windows |
| Deadspace | < 6 cc                                     |

### CO<sub>2</sub> DISPLAY

|                    |   |
|--------------------|---|
| Screen display     | CO <sub>2</sub> waveform and ETCO <sub>2</sub> and INCO <sub>2</sub> numerics                                 |
| Measurement ranges | ETCO <sub>2</sub> : 0-99 mmHg, 0-13 kPa, 0-23%<br>INCO <sub>2</sub> : 0-25 mmHg, 0-5 kPa, 0-5%                |
| Display ranges     | ETCO <sub>2</sub> and INCO <sub>2</sub> same as measurement range   |
| Units              | mmHg, kPa, %; user-selectable   |
| Sweep speed        | 3.13, 6.25, 12.5 mm/sec; user-selectable  |
| Response modes     | FAST: 15 sec sampling time period<br>NORMAL: 30 sec sampling time period<br>SLOW: 45 sec sampling time period |
| Gas compensation   | see User's Guide  |
| Alarm limit ranges | ETCO <sub>2</sub> : 0-99 mmHg, 0-14 kPa, 0-14%<br>INCO <sub>2</sub> : 2-25 mmHg                               |
| Resolution         | 1 mmHg  |

|   |  |
|---|--|
| Accuracy                                | +/- 3 mmHg (0-30 mmHg CO <sub>2</sub> )<br>+/- 10% of reading (31-99 mmHg CO <sub>2</sub> )  |
| Waveform rise time                      | 120 ms maximum   |
| Altitude error                          | +/- .4%/1000 ft  |
| <b>BREATH RATE DISPLAY</b>              |  |
| Screen display                          | numeric  |
| Units                                   | bpm  |
| Range                                   | 2-150 bpm  |
| <b>APNEA ALARMS AND TICKETS</b>         |  |
| Apnea ticket                            | set to auto print after apnea event and after 1 minute continued apnea   |
| Apnea alarm accuracy                    | +/- 1 sec  |
| Resolution                              | 5 sec  |
| Alarm limits range, adult and pediatric | 15-30 sec delay, 5 sec increments  |
| <b>BAROMETRIC PRESSURE</b>              |  |
| Pressure compensation                   | automatic  |
| Operating range                         | -2000 ft to 15,000 ft  |
| Screen display                          | numeric (CO <sub>2</sub> status window)  |
| Units                                   | mmHg or kPa  |
| Accuracy                                | +/- 3mmHg  |
| <b>IN-SERVICE VALUES</b>                |  |
| ETCO <sub>2</sub>                       | initial value: 38, alternate value: 60   |
| INCO <sub>2</sub>                       | initial value: 0, alternate value: 8   |
| Breath rate                             | initial value: 12, alternate value: 31   |
| <b>TEMPERATURE</b>                      |  |
| RANGE                                   | 17°C to 50°C; 62.6°F to 122°F  |
| DISPLAYS                                | T1   |
| PROBES                                  | Compatible with YSI Series 400 and 700 and Electromedics Series 2100 probes. HP side panel only compatible with YSI 400 and has HP connector |
| UNITS                                   | °C or °F, user selectable  |
| ACCURACY                                | +/- .1°C (+/- .2°F) plus probe tolerance   |
| RESOLUTION                              | .1°C or °F   |
| ELECTROSURGERY SUPPRESSION              | Included in all EL display monitors  |
| <b>ALARMS</b>                           |  |
| INDICATORS                              | ALARM light, ALARM(S) OFF light, Audible tone<br>Lights continually flash 0.5 sec on and 0.5 sec off if an alarm is suspected                |
| TONE FREQUENCY                          | 900 Hz<br>Tone is steady for a patient alarm and sounds for 1- sec every 4 sec for an equipment alert  |

SELECTABLE TONE VOLUME  
LIMITS  
CONTROL  
ALARM ON TACHYCARDIAS

low, medium, high  
settable on all parameters  
Automatic preset or manual settings  
Most tachycardias will alarm in less than 8  
sec.

## DISPLAY

### GENERAL

Matrix  
Active viewing area

552 X 256 pixels  
145.75 mm X 67.56 mm

### ELECTROLUMINESCENT DISPLAY

Viewing angle  
Display window  
Display color  
Display background color

> 160° Horizontal and vertical  
contrast enhancement filter  
amber  
black

## MONITOR (Environmental)

OPERATING TEMPERATURE  
SHIPPING AND STORAGE TEMPERATURE  
OPERATING ALTITUDE  
SHIPPING AND STORAGE ALTITUDE  
OPERATING RELATIVE HUMIDITY  
SHIPPING AND STORAGE RELATIVE  
HUMIDITY  
SHOCK  
VIBRATION

0-40°C  
-20-60°C  
-2000-15000 ft  
-2000-40000 ft  
15-95%, noncondensing  
15-95%, noncondensing  
50 g  
Random vibration, .02 g<sup>2</sup>/Hz from 10 - 500 Hz,  
ramping down to .002 g<sup>2</sup>/Hz at 2000 Hz.  
Operating 1 hr per axis, 3 hr per test.  
Per IEC 601-1-2

### ELECTROMAGNETIC INTERFERENCE

## MONITOR (Physical)

### PROTECTION CLASSIFICATIONS

Type of protection against electric shock  
- monitor powered by power adapter  
Degree of protection against electric shock  
Method of disinfection  
Flammable anesthetics

- Class I (protectively earthed)  
- Type CF equipment, Defibrillator-proof  
- not suitable for autoclaving  
- not suitable for use with flammable anesthetics

### MONITOR ONLY

Height  
Width  
Depth  
Weight

6.65 in  
8.25 in  
5.10 in  
6.25 lb

## MONITOR WITH EXPANSION MODULE

|        |          |
|--------|----------|
| Height | 9.65 in  |
| Width  | 8.25 in  |
| Depth  | 7.50 in  |
| Weight | 12.68 lb |

## PRINTER

### OPERATION

|                      |   |
|----------------------|---|
| Operating modes      | Continuous, Snapshot, Freeze Print, Auto Interval Print, Auto Interval Trend, Tabular Trend, Alarm Print, Cuff Ticket, Apnea Ticket                                     |
| Auto Print Intervals | 15 min, 30 min, 1 hr, 2 hr, 4 hr  |
| Auto trend shifts    | once every 4 hr   |
| Number of waveforms  | up to three: ECG, SpO <sub>2</sub> , P1, P2, CO <sub>2</sub>  |
| Grid                 | 5 mm and 1 mm gradations  |
| Annotation           | Date, Time, Print mode, Speed, Heart rate, Systolic, Diastolic, Mean, SpO <sub>2</sub> , Breath rate, ETCO <sub>2</sub> , INCO <sub>2</sub> , Temperature, Pacer status |
| Printing Speeds      | 6.25, 12.5, 25.0 mm/sec, simulated 6.25   |

### PRINTER MECHANISM

|                      |                                |
|----------------------|--------------------------------|
| Printing method      | thermally sensitive dot method |
| Dot structure        | 320 dots per line              |
| Printing width       | 53 mm                          |
| Horizontal dot pitch | 0.165 mm, 6 dots/mm            |
| Vertical dot pitch   | 0.165 mm                       |
| Paper feed method    | friction feed                  |
| Paper feed precision | +/- 2% @ 25°C and 60% RH       |
| Paper width          | 60 mm                          |
| Reliability          | 30 million pulses/dot          |

### ENVIRONMENTAL

|                                     |   |
|-------------------------------------|---|
| Operating temperature               | 5-40°C  |
| Shipping and storage temperature    | -20-60°C  |
| Operating relative humidity         | 35% to 85% noncondensing  |
| Shipping, storage relative humidity | 15% to 90% noncondensing  |
| Shipping and operating altitude     | -2000 to 15000 ft   |
| Storage altitude                    | -2000 to 40000 ft   |
| Shock                               | 30 g  |
| Vibration                           | Random vibration, .02 g <sup>2</sup> /Hz from 10 to 500 Hz, ramping down to .002 g <sup>2</sup> /Hz at 2000 Hz. Operating 1 hour per axis, 3 hours per test. per standard IEC 601-1-2 |
| EMI                                 |   |

### PAPER

|  |                                   |
|--|-----------------------------------|
| Short-term storage environment<br>(up to 7 days) | -20-40°C, 5% to 80% noncondensing |
| Long-term storage environment<br>(up to 5 years) | 25°C (optimal), 65% noncondensing |



## POWER

|                                       |   |
|---------------------------------------|---|
| MODE OF OPERATION                     | Continuous  |
| BATTERY PACK TYPE                     | Sealed gel-type lead acid                                       |
| BATTERY PACK CAPACITY                 | Monitor only - 8 volts, 3 amp-hours                             |
|                                       | Monitor with expansion modules - 8 volts, 6 amp-hours           |
| BATTERY RECHARGER CIRCUITRY           | Internal, powered by external power adapter                     |
| DC INPUT POWER REQUIRED               | 12-28 Volts, 10.5 Watts, w/CO <sub>2</sub> : 25 Watts           |
| INPUT FUSE RATING                     | 3 A/250V, Slow-blow, Type 2AG (.57 X .177 in.)                  |
| OPERATING TIMES ON BATTERY            | Range of 3.5 to 4.5 hr depending on product configuration       |
| BATTERY RECHARGE TIME WITH 206 EL ON  | Range of 8 to 12 hr typical, depending on product configuration |
| BATTERY RECHARGE TIME WITH 206 EL OFF | Range of 6 to 8 hr depending on product configuration           |

## POWER ADAPTERS

### UNIVERSAL POWER ADAPTER, PART NO.

503-0054-00

|                           |  |
|---------------------------|--|
| Length                    | 5.0 in.  |
| Width                     | 3.6 in.  |
| Height                    | 3.1 in.  |
| Weight                    | 3.1 lb   |
| Rated input               | 100 V-120 VAC, 500 mA, 50/60 Hz                  |
| Rated fuses               | T800 mA/250 V, Time-delay, 5 X 20mm              |
| Rated output (continuous) | 16-24 VDC, 25 VA                                 |
| Additional Features       | Detachable power cord, pilot light, mains switch |